


Genotyping Report

| | | | | | |
|-----------|-----------|-------------|-----------|--------------------|-------------|
| Strain ID | T005371 | Strain Type | CKO(Cas9) | Genetic Background | C57BL/6JGpt |
| Designer | Ya'nan Xu | Gene Name | Wwp2 | | |

1. Strategy of Genotyping

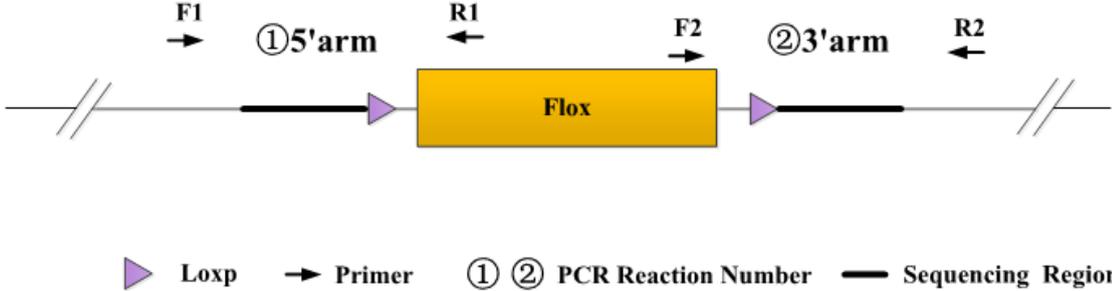


Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.
 Heterozygote: ①PCR reaction obtains a WT band and a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.
 Homozygote: ①PCR reaction obtains a single Targeted band; ②PCR reaction obtains a single Targeted band.
 Note: The sizes of WT and Targeted band are shown below.

2. Primer Information

| PCR No. | Primer No. | Sequence | Band Size |
|----------|------------|---------------------------|-----------------------------|
| ①(5'arm) | T005371-F1 | AGCACTCCTGTCTTCCTCAGCT | WT: 264bp Targeted:369bp |
| | T005371-R1 | CACTGCCTCTGAGCATCTCCAT | |
| ②(3'arm) | T005371-F2 | GGAGGAAGAGGATGATGTGTTTCAT | WT: 280bp Targeted:386bp |
| | T005371-R2 | GCAGTCAATAAGCTCTGGTCACCT | |

3. Gel Image & Conclusion

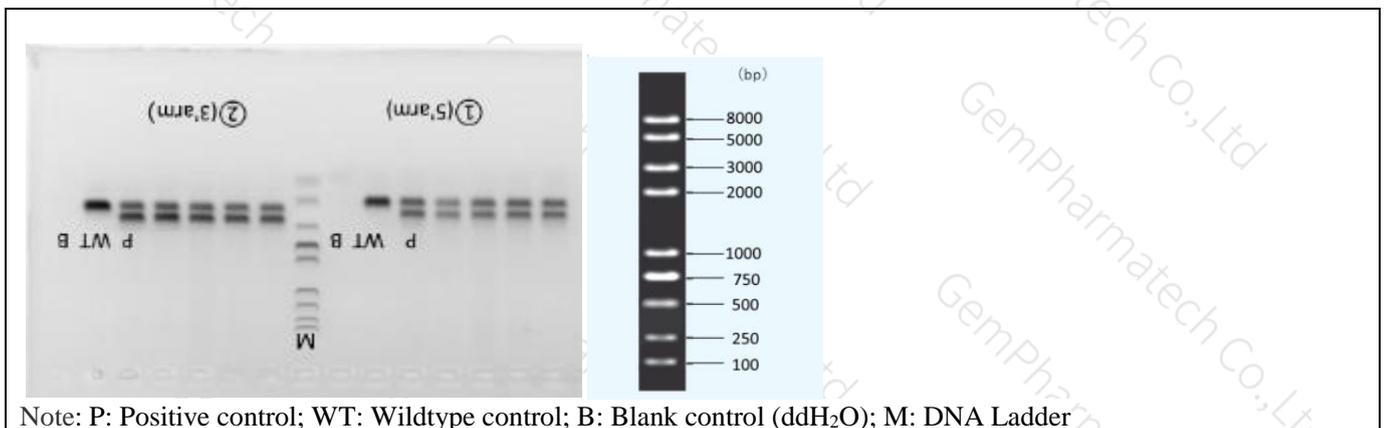

- ① Control (WT) : It is an important reference mark for whether the PCR reaction is successful and whether the product band position and size meet the theoretical requirements.
- ② Control (B) : PCR amplification was performed without template in the PCR reagent to monitor whether the reagent was contaminated.

4. PCR Condition

| PCR Reaction Component | | | |
|----------------------------------|--|------|-------------|
| Seg. | reaction component | | Volume (μl) |
| 1 | 2 × Rapid Taq Master Mix (Vazyme P222) | | 12.5 |
| 2 | ddH ₂ O | | 9.5 |
| 3 | Primer A(10pmol/μl) | | 1 |
| 4 | Primer B(10pmol/μl) | | 1 |
| 5 | Template(≈100ng/μl) | | 1 |
| PCR program ① priority selection | | | |
| Seg. | Temp. | Time | Cycle |
| 1 | 95°C | 5min | |
| 2 | 98°C | 30s | 20× |
| 3 | 65°C* (-0.5°C/cycle) | 30s | |
| 4 | 72°C | 45s* | |
| 5 | 98°C | 30s | 20× |
| 6 | 55°C* | 30s | |
| 7 | 72°C | 45s* | |
| 8 | 72°C | 5min | |
| 9 | 10°C | hold | |
| PCR program ② the second choice | | | |
| Seg. | Temp. | Time | Cycle |
| 1 | 95°C | 5min | |
| 2 | 98°C | 30s | 35× |
| 3 | 58°C* | 30s | |
| 4 | 72°C | 45s* | |
| 5 | 72°C | 5min | |
| 6 | 10°C | hold | |

Note*: Annealing temperature and extension time can be determined according to the actual amplification situation and amplification enzyme efficiency.

