


江苏集萃药康生物科技股份有限公司 GemPharmatech Co.,Ltd

		Genoty	oing Report		- Kr
Strain ID	T027129	Strain Type	CKO(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Zifan Lin	Gene Name	~~~~~~~~~~~~~~~~~~~~~~~~~~~~~~~~~~~~~~	Dcaf10	6
94	35	10	· · · · · · · · · · · · · · · · · · ·	Qr.	3 Cx

1. Strategy of Genotyping

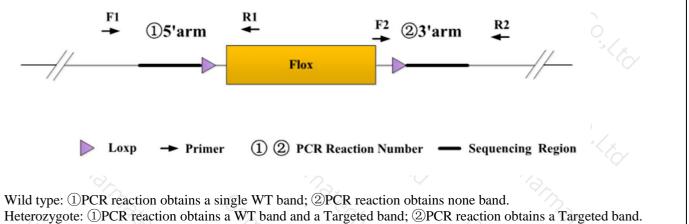


Heterozygote: (1)PCR reaction obtains a W1 band and a Targeted band; (2)PCR reaction obtains a Targeted band, Homozygote: (1)PCR reaction obtains a single Targeted band; (2)PCR reaction obtains a Targeted band. Note: The sizes of WT and Targeted band are shown below.

2. Primer Information

PCR No.	Primer No.	Sequence	Band Size
(1)(5'arm)	T027129-F1 GGCTACACAGAGAAACCCTGTCTCA A		WT: 324 bp Targeted: 429 bp
	T027129-R1	129-R1 CCAAGCACTCTATAATTTGGGGTGC	
(2)(3'arm)	T027129-F2	T027129-F2 ATCGCATTGTCTGAGTACGTG	
	T027129-R2	TAGCCCTTCTTCTCTTTCAGGTTTCC	Targeted: 260 bp

3. Gel Image & Conclusion

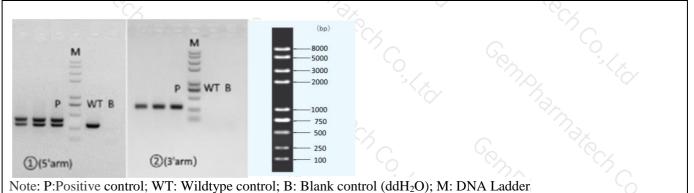

① Control (WT) : It is an important reference mark for whether the PCR reaction is successful and whether the



江苏集萃药康生物科技股份有限公司 GemPharmatech Co.,Ltd

product band position and size meet the theoretical requirements.

⁽²⁾ Control (B) : PCR amplification was performed without template in the PCR reagent to monitor whether the reagent was contaminated.

4. PCR Condition

PCR Reaction	Component))	AN CA		
Seg.	reaction	reaction component			
1 75,	2 × Rapid Taq Master Mix (Vazyr	2 × Rapid Taq Master Mix (Vazyme P222)			
2	ddH2O	$\sim \sim \sim$	9.5		
3	Primer A(10pmol/µl)		1 A		
4	Primer B(10pmol/μl)	Primer B(10pmol/µl)			
5	Template(≈100ng/µl)	Template(≈100ng/μl)			
PCR program	① priority selection				
Seg.	Temp.	Time	Cycle		
1	95°C	5min	CIAFIN -		
2 6	98°C	30s	20×		
3	65℃*(-0.5℃/cycle)	30s	$\overline{\mathcal{D}}_{\mathcal{A}}$		
4	72℃	45s*			
5	98°C	30s	20×		
6	55℃*	30s			
7 2	72°C	45s*			
8	72°C	5min	Par ila		
9	10 °C	hold	173× 19		
PCR program	② the second choice	- 73× 60	202 202		
Seg.	Temp.	Time	Cycle		
1	95°C	5min			
2	98°C	30s	35×		
3	58°C*	58℃* 30s			
4	72°C	45s*			
5	72℃	5min			
6	10°C	hold	ng p		

Note*: Annealing temperature and extension time can be determined according to the actual amplification situation and amplification enzyme efficiency.