


## Genotyping Report

|           |                 |             |              |                    |             |
|-----------|-----------------|-------------|--------------|--------------------|-------------|
| Strain ID | JS25058-T024764 | Strain Type | CKO(Cas9)    | Genetic Background | C57BL/6JGpt |
| Designer  | Ya'nan Xu       | Gene Name   | <i>Pcmd1</i> |                    |             |

### 1. Strategy of Genotyping

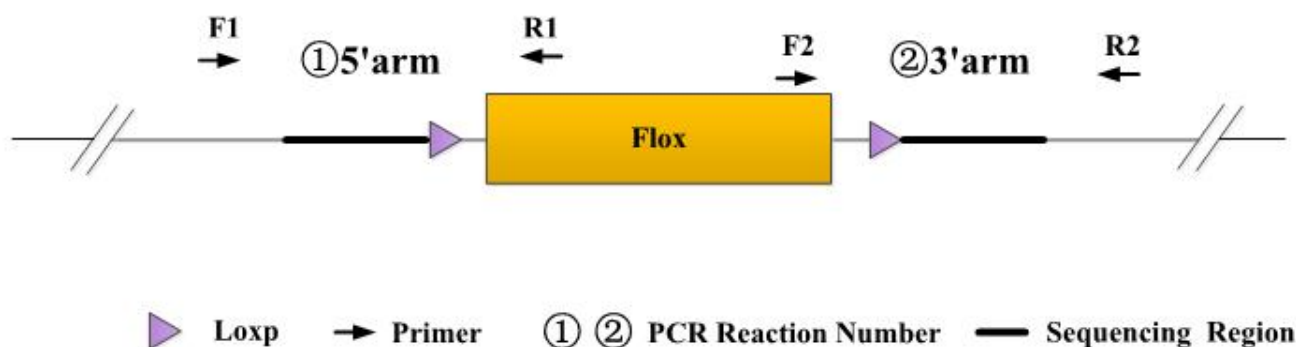

Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.

Heterozygote: ①PCR reaction obtains a WT band and a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.

Homozygote: ①PCR reaction obtains a single Targeted band; ②PCR reaction obtains a single Targeted band.

Note: The sizes of WT and Targeted band are shown below.

### 2. Primer Information

| PCR No.  | Primer No.     | Sequence                   | Band Size                       |
|----------|----------------|----------------------------|---------------------------------|
| ①(5'arm) | T024764(P1)-F1 | AGGTCTTTCAGCCATCATCTTG     | WT: 274bp<br>Targeted:<br>379bp |
|          | T024764(P1)-R1 | CCTCCAGCCCTCTGTTTACTAACTTG |                                 |
| ②(3'arm) | T024764(P1)-F2 | GCGATCCCTGCTATACAGTAAGAAAC | WT: 291bp<br>Targeted:<br>397bp |
|          | T024764(P1)-R2 | CAACCCACATGACACATACAAATG   |                                 |

### 3. Gel Image & Conclusion



Note: P: Heterozygous samples; WT: Wildtype control; B: Blank control (ddH<sub>2</sub>O); M: DNA Ladder

- ① Control (WT) : It is an important reference mark for whether the PCR reaction is successful and whether the product band position and size meet the theoretical requirements.
- ② Control (B) : PCR amplification was performed without template in the PCR reagent to monitor whether the reagent was contaminated.

#### 4. PCR Condition

| PCR Reaction Component           |                                        |      |             |
|----------------------------------|----------------------------------------|------|-------------|
| Seg.                             | reaction component                     |      | Volume (μl) |
| 1                                | 2 × Rapid Taq Master Mix (Vazyme P222) |      | 12.5        |
| 2                                | ddH2O                                  |      | 9.5         |
| 3                                | Primer A(10pmol/μl)                    |      | 1           |
| 4                                | Primer B(10pmol/μl)                    |      | 1           |
| 5                                | Template(20~80ng/μl)                   |      | 1           |
| PCR program I priority selection |                                        |      |             |
| Seg.                             | Temp.                                  | Time | Cycle       |
| 1                                | 95℃                                    | 5min | 20×         |
| 2                                | 98℃                                    | 30s  |             |
| 3                                | 65℃* (-0.5℃/cycle)                     | 30s  |             |
| 4                                | 72℃                                    | 45s* |             |
| 5                                | 98℃                                    | 30s  | 15×         |
| 6                                | 55℃*                                   | 30s  |             |
| 7                                | 72℃                                    | 45s* |             |
| 8                                | 72℃                                    | 5min |             |
| 9                                | 10℃                                    | hold |             |
| PCR program II the second choice |                                        |      |             |
| Seg.                             | Temp.                                  | Time | Cycle       |
| 1                                | 95℃                                    | 5min | 35×         |
| 2                                | 98℃                                    | 30s  |             |
| 3                                | 58℃*                                   | 30s  |             |
| 4                                | 72℃                                    | 45s* |             |
| 5                                | 72℃                                    | 5min |             |
| 6                                | 10℃                                    | hold |             |

Note\*: Annealing temperature and extension time can be determined according to the actual amplification situation and amplification enzyme efficiency.