


## Genotyping Report

|           |           |             |              |                    |             |
|-----------|-----------|-------------|--------------|--------------------|-------------|
| Strain ID | T051789   | Strain Type | CKO(Cas9)    | Genetic Background | C57BL/6JGpt |
| Designer  | Ya'nan Xu | Gene Name   | <i>Senp5</i> |                    |             |

### 1. Strategy of Genotyping

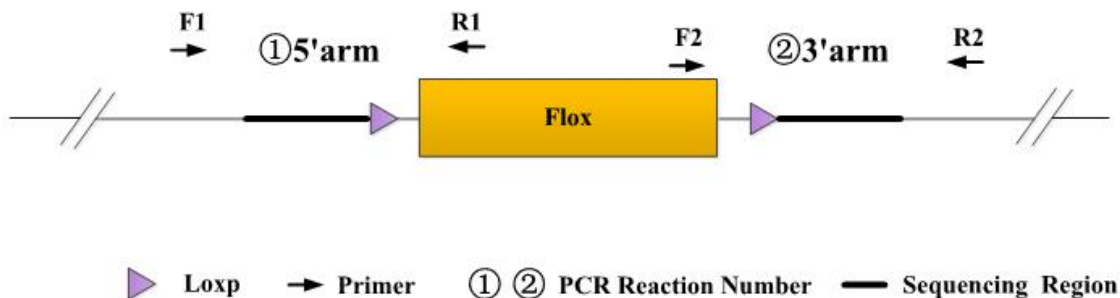

Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.

Heterozygote: ①PCR reaction obtains a WT band and a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.

Homozygote: ①PCR reaction obtains a single Targeted band; ②PCR reaction obtains a single Targeted band.


Note: The sizes of WT and Targeted band are shown below.

### 2. Primer Information

| PCR No.  | Primer No. | Primer Name    | Sequence                    | Band Size                   |
|----------|------------|----------------|-----------------------------|-----------------------------|
| ①(5'arm) | F1         | T051789(P3)-F1 | AGATGGTTGTGAGCCACCATGT      | WT: 332bp<br>Targeted:437bp |
|          | R1         | T051789(P3)-R1 | TCCTACACTCTGGTGTTACAGGTGTGT |                             |
| ②(3'arm) | F2         | T051789(P3)-F2 | ATAACAGCAGATTGCTGTGGCCTA    | WT: 243bp<br>Targeted:349bp |
|          | R2         | T051789(P3)-R2 | TTCAGTGCTAGAGACTGTACACAGGCT |                             |

### 3. Gel Image & Conclusion

|  |
|--|
|  |
|--|

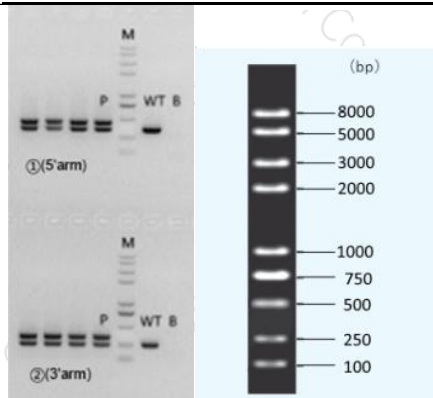

Note: P: Positive control; WT: Wildtype control; B: Blank control (ddH<sub>2</sub>O); M: DNA Ladder

① Control (WT) : It is an important reference mark for whether the PCR reaction is successful and whether the product band position and size meet the theoretical requirements.

② Control (B) : PCR amplification was performed without template in the PCR reagent to monitor whether the reagent was contaminated.

#### 4. PCR Condition

| PCR Reaction Component           |                                        |      |             |
|----------------------------------|----------------------------------------|------|-------------|
| Seg.                             | reaction component                     |      | Volume (μl) |
| 1                                | 2 × Rapid Taq Master Mix (Vazyme P222) |      | 12.5        |
| 2                                | ddH2O                                  |      | 9.5         |
| 3                                | Primer A(10pmol/μl)                    |      | 1           |
| 4                                | Primer B(10pmol/μl)                    |      | 1           |
| 5                                | Template(≈100ng/μl)                    |      | 1           |
| PCR program ① priority selection |                                        |      |             |
| Seg.                             | Temp.                                  | Time | Cycle       |
| 1                                | 95℃                                    | 5min | 20×         |
| 2                                | 98℃                                    | 30s  |             |
| 3                                | 65℃* ( -0.5℃/cycle )                   | 30s  |             |
| 4                                | 72℃                                    | 45s* |             |
| 5                                | 98℃                                    | 30s  | 20×         |
| 6                                | 55℃*                                   | 30s  |             |
| 7                                | 72℃                                    | 45s* |             |
| 8                                | 72℃                                    | 5min |             |
| 9                                | 10℃                                    | hold |             |
| PCR program ② the second choice  |                                        |      |             |
| Seg.                             | Temp.                                  | Time | Cycle       |

|   |       |      |     |
|---|-------|------|-----|
| 1 | 95℃   | 5min |     |
| 2 | 98℃   | 30s  | 35× |
| 3 | 58℃ * | 30s  |     |
| 4 | 72℃   | 45s* |     |
| 5 | 72℃   | 5min |     |
| 6 | 10℃   | hold |     |

Note\*: Annealing temperature and extension time can be determined according to the actual amplification situation and amplification enzyme efficiency.