




## Genotyping Report

Strain ID	T022938	Strain Type	CKO(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Ya'nan Xu	Gene Name			<i>Ndufb6</i>

### 1. Strategy of Genotyping

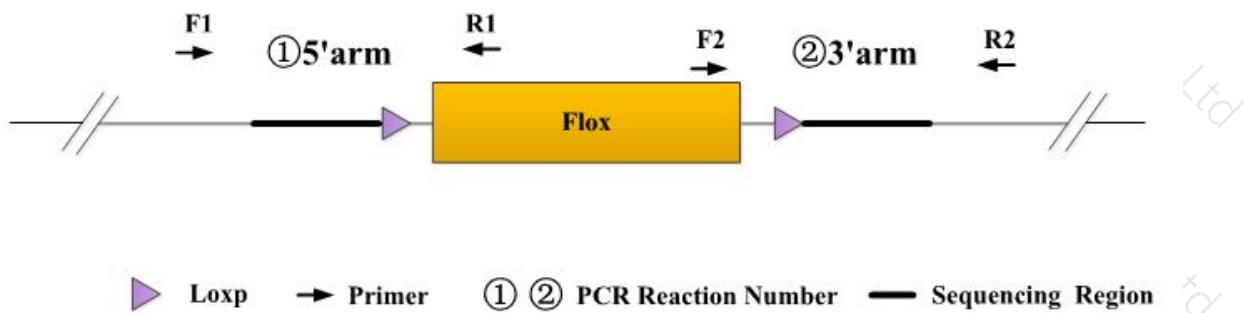

Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.

Heterozygote: ①PCR reaction obtains a WT band and a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.

Homozygote: ①PCR reaction obtains a single Targeted band; ②PCR reaction obtains a single Targeted band.


Note: The sizes of WT and Targeted band are shown below.

### 2. Primer Information

PCR No.	Primer No.	Primer Name	Sequence	Band Size
①(5'arm)	F1	JS23232-Ndufb6-5wt-tF1	GGATGTGTTCTAGGTAGCAGAAGC	WT:368bp Targeted:473bp
	R1	JS23232-Ndufb6-5wt-tR1	AATTCTGATGCTGTGGCTGCTC	
②(3'arm) GC:36.8%	F2	JS23232-Ndufb6-3wt-tF1A	CTGAGGTTCTGTCATTCCCAGATT	WT:389bp Targeted:495bp
	R2	JS23232-Ndufb6-3wt-tR1A	CTGGGGTCTGACATACACTGATTACA	

### 3. Gel Image & Conclusion



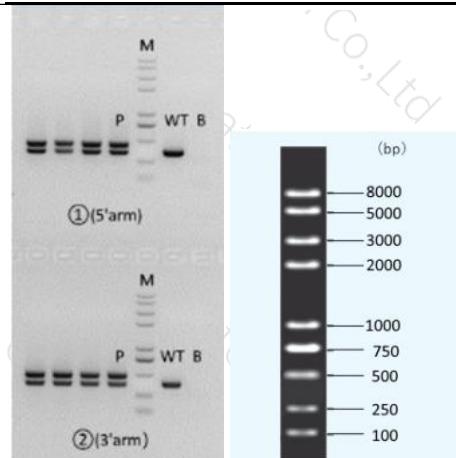

(bp)  
8000  
5000  
3000  
2000  
1000  
750  
500  
250  
100

Note: P: Heterozygous samples; WT: Wildtype control; B: Blank control (ddH<sub>2</sub>O); M: DNA Ladder

① Control (WT) : It is an important reference mark for whether the PCR reaction is successful and whether the product band position and size meet the theoretical requirements.

② Control (B) : PCR amplification was performed without template in the PCR reagent to monitor whether the reagent was contaminated.

#### 4. PCR Condition

PCR Reaction Component			
Seg.	reaction component	Volume (μl)	
1	2 × Rapid Taq Master Mix (Vazyme P222)	12.5	
2	ddH <sub>2</sub> O	9.5	
3	Primer A(10pmol/μl)	1	
4	Primer B(10pmol/μl)	1	
5	Template(20~80ng/μl)	1	

PCR program I priority selection			
Seg.	Temp.	Time	Cycle
1	95 °C	5min	
2	98 °C	30s	20×
3	65 °C * (-0.5 °C/cycle)	30s	
4	72 °C	45s*	15×
5	98 °C	30s	
6	55 °C *	30s	
7	72 °C	45s*	
8	72 °C	5min	
9	10 °C	hold	

PCR program II the second choice			
Seg.	Temp.	Time	Cycle



江苏集萃药康生物科技股份有限公司  
GemPharmatech Co.,Ltd

1	95°C	5min	
2	98°C	30s	35x
3	58°C *	30s	
4	72°C	45s*	
5	72°C	5min	
6	10°C	hold	

Note\*: Annealing temperature and extension time can be determined according to the actual amplification situation and amplification enzyme efficiency.