


Genotyping Report

Strain ID	T016217	Strain Type	CKO(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Tiantian Sun	Gene Name	<i>Rnfl25</i>		

1. Strategy of Genotyping

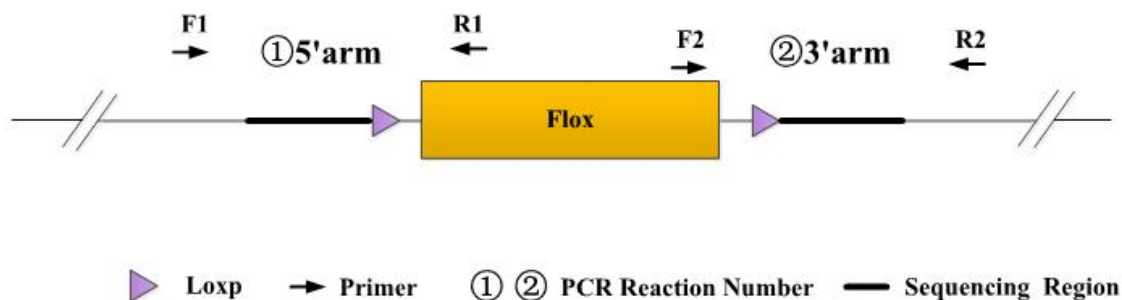

Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.

Heterozygote: ①PCR reaction obtains a WT band and a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.


Homozygote: ①PCR reaction obtains a single Targeted band; ②PCR reaction obtains a single Targeted band.

Note: The sizes of WT and Targeted band are shown below.

2. Primer Information

PCR No.	Primer No.	Sequence	Band Size
①(5'arm)	T016217(P3)-F1	CTTTGCCTCAGAGCCTGGAAA	WT:267bp Targeted:372bp
	T016217(P3)-R1	CCTGTGTGTCCTTGCCTCAAA	
②(3'arm) GC% 35.7	T016217(P3)-F2	CCTAGAGATCCAAAGGATCACTCTGG	WT:362bp Targeted:471bp
	T016217(P3)-R2	TCCTTCTACCTTGGTTTCCTGCA	

3. Gel Image & Conclusion

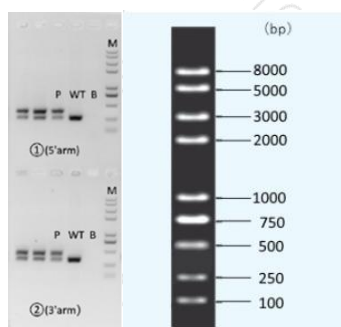

Note: P: Positive control; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

- ① Control (WT) : It is an important reference mark for whether the PCR reaction is successful and whether the product band position and size meet the theoretical requirements.
- ② Control (B) : PCR amplification was performed without template in the PCR reagent to monitor whether the reagent was contaminated.

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Rapid Taq Master Mix (Vazyme P322)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program ① priority selection			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃* (-0.5℃/cycle)	30s	
4	72℃	45s*	20×
5	98℃	30s	
6	55℃*	30s	
7	72℃	45s*	
8	72℃	5min	
9	10℃	hold	
PCR program ② the second choice			
Seg.	Temp.	Time	Cycle
1	95℃	5min	35×
2	98℃	30s	
3	58℃*	30s	
4	72℃	45s*	
5	72℃	5min	
6	10℃	hold	

Note*: Annealing temperature and extension time can be determined according to the actual amplification situation and amplification enzyme efficiency.



江苏集萃药康生物科技股份有限公司
GemPharmatech Co.,Ltd