


Genotyping Report

| | | | | | |
|-----------|--------------|-------------|---------------|--------------------|-------------|
| Strain ID | T005805 | Strain Type | CKO(Cas9) | Genetic Background | C57BL/6JGpt |
| Designer | Tiantian Sun | Gene Name | <i>Gpr137</i> | | |

1. Strategy of Genotyping

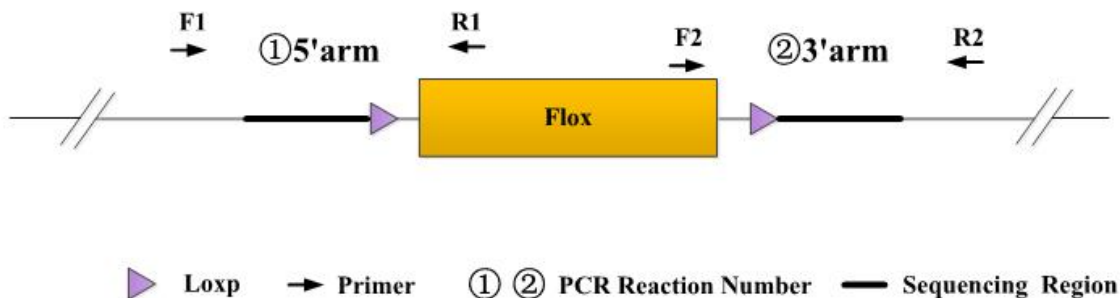

Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.

Heterozygote: ①PCR reaction obtains a WT band and a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.

Homozygote: ①PCR reaction obtains a single Targeted band; ②PCR reaction obtains a single Targeted band.

Note: The sizes of WT and Targeted band are shown below.

2. Primer Information

| PCR No. | Primer No. | Primer Name | Sequence | Band Size |
|----------|------------|-------------|-------------------------|----------------------------|
| ①(5'arm) | F1 | T005805-F1 | CACCTGTGGTGGAGAAAGATGG | WT:285bp Targeted:366bp |
| | R1 | T005805-R1 | GGGTGATTGAGAGCTACAGCATC | |
| ②(3'arm) | F2 | T005805-F2 | GAAGCGAGAGGTTTGTCTGGT | WT:474bp Targeted:552bp |
| | R2 | T005805-R2 | GGACAGAAGGAGATAATTGCCTC | |

3. Gel Image & Conclusion



Note: P: Heterozygous samples; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

① Control (WT) : It is an important reference mark for whether the PCR reaction is successful and whether the

product band position and size meet the theoretical requirements.

② Control (B) : PCR amplification was performed without template in the PCR reagent to monitor whether the reagent was contaminated.

4. PCR Condition

Generally recommend to use Vazyme P222; If the sequences contain special structures such as GC% $\geq 60\%$ or GC% $\leq 40\%$, recommend to use Vazyme P515

| PCR Reaction Component | | | |
|----------------------------------|--|------|-------------|
| Seg. | reaction component | | Volume (μl) |
| 1 | 2 × Rapid Taq Master Mix(Vazyme P222) or 2 × Phanta Max Master Mix (Vazyme P515) | | 12.5 |
| 2 | ddH2O | | 9.5 |
| 3 | Primer A(10pmol/μl) | | 1 |
| 4 | Primer B(10pmol/μl) | | 1 |
| 5 | Template(20~80ng/μl) | | 1 |
| PCR program I priority selection | | | |
| Seg. | Temp. | Time | Cycle |
| 1 | 95℃ | 5min | 20× |
| 2 | 98℃ | 30s | |
| 3 | 65℃*（-0.5℃/cycle） | 30s | |
| 4 | 72℃ | 45s* | |
| 5 | 98℃ | 30s | 15× |
| 6 | 55℃* | 30s | |
| 7 | 72℃ | 45s* | |
| 8 | 72℃ | 5min | |
| 9 | 10℃ | hold | |
| PCR program II the second choice | | | |
| Seg. | Temp. | Time | Cycle |
| 1 | 95℃ | 5min | 35× |
| 2 | 98℃ | 30s | |
| 3 | 58℃* | 30s | |
| 4 | 72℃ | 45s* | |
| 5 | 72℃ | 5min | |
| 6 | 10℃ | hold | |

Note*: Annealing temperature and extension time can be determined according to the actual amplification situation and amplification enzyme efficiency.