

Genotyping Report

Strain ID	T009208	Strain Type	CKO(Cas9)	Genetic Background	C57BL/6J
Designer	Chen Chen	Gene Name	Cul4b		

1. Strategy of Genotyping



Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.

Heterozygote: ①PCR reaction obtains a WT band and a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.

Homozygote: ①PCR reaction obtains a single Targeted band; ②PCR reaction obtains a single Targeted band.


Note: The sizes of WT and Targeted band are shown below.

2. Primer Information

PCR No.	Primer No.	Primer Name	Sequence	Band Size
①(5'arm) GC:40%	F1	JS03084-Cul4b-5wt-tF1	GCATAGGTTCTGGAGTCATTTAGGTG	WT:351bp Targeted:453bp
	R1	JS03084-Cul4b-5wt-tR1	CAGAAGCTGACATGAATCTTGTAGGG	
②(3'arm) GC:39.4%	F2	JS03084-Cul4b-3wt-tF1	CCTCATTTGGTTCTTTGGTCCCAG	WT:273bp Targeted:376bp
	R2	JS03084-Cul4b-3wt-tR1	CGTAGTTTGACCTTCCCTATTCAACC	

3. Gel Image & Conclusion

--

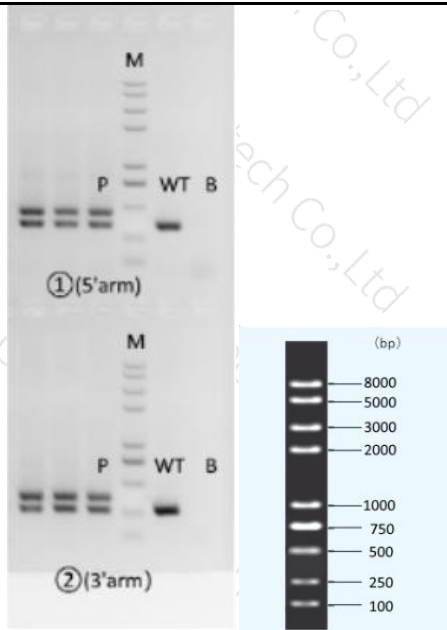

Note: P: Heterozygous samples; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

① Control (WT) : It is an important reference mark for whether the PCR reaction is successful and whether the product band position and size meet the theoretical requirements.

② Control (B) : PCR amplification was performed without template in the PCR reagent to monitor whether the reagent was contaminated.

4. PCR Condition

(Generally recommend to use Vazyme P222; If the sequences contain special structures such as GC% $\geq$ 60% or GC% $\leq$ 40%, recommend to use Vazyme P515.)

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Rapid Taq Master Mix(Vazyme P222) or 2 × Phanta Max Master Mix (Vazyme P515)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(20~80ng/μl)		1
PCR program I priority selection			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃* (-0.5℃/cycle)	30s	
4	72℃	45s*	
5	98℃	30s	15×

6	55℃ *	30s	
7	72℃	45s*	
8	72℃	5min	
9	10℃	hold	
PCR program II the second choice			
Seg.	Temp.	Time	Cycle
1	95℃	5min	
2	98℃	30s	35×
3	58℃ *	30s	
4	72℃	45s*	
5	72℃	5min	
6	10℃	hold	

Note*: Annealing temperature and extension time can be determined according to the actual amplification situation and amplification enzyme efficiency.