




Genotyping Report

Strain ID	T017865	Strain Type	CKO(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Ya'nan Xu	Gene Name			<i>Jaml</i>

1. Strategy of Genotyping

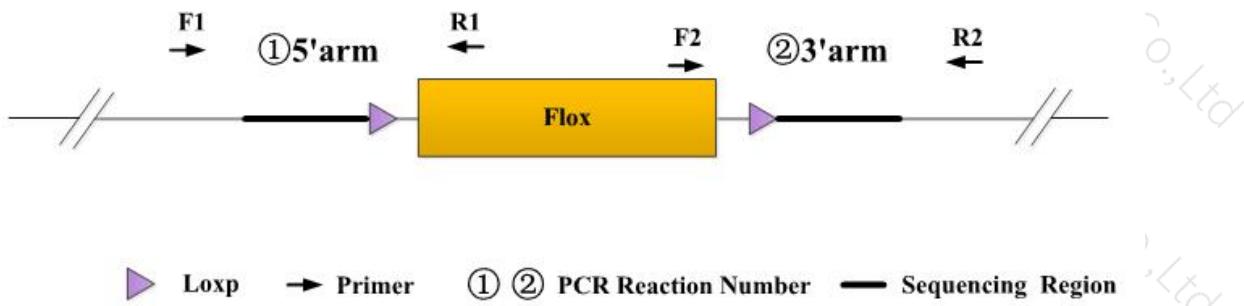

Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.

Heterozygote: ①PCR reaction obtains a WT band and a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.

Homozygote: ①PCR reaction obtains a single Targeted band; ②PCR reaction obtains a single Targeted band.


Note: The sizes of WT and Targeted band are shown below.

2. Primer Information

PCR No.	Primer No.	Sequence	Band Size
①(5'arm) Alternative	T017865-F1A	CCGATGTTGTAAGTCAGTGACAGG	WT:402bp Targeted:507bp
	T017865-R1A	AGATACAACACATTCTCTAACCTGGGT	
②(3'arm) Alternative	T017865-F2A	CAATCTCATTCTTAGCCAAGTGTGG	WT:416bp Targeted:517bp
	T017865-R2A	GGACCCTTCTGTCTAACCTGTAC	

3. Gel Image & Conclusion




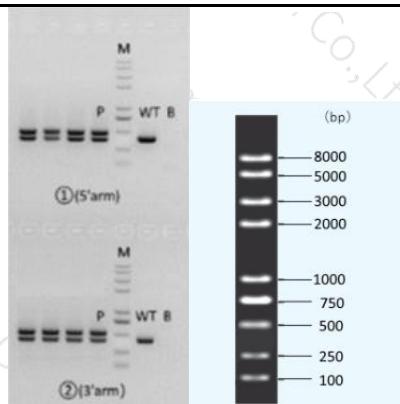

GemPharmatech Co.,Ltd

GemPharmatech Co.,Ltd

GemPharmatech Co.,Ltd

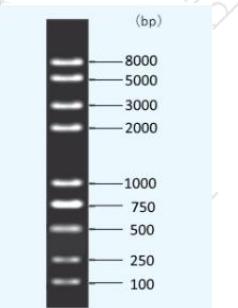

GemPharmatech Co.,Ltd

GemPharmatech Co.,Ltd

Note: P: Positive control; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

① Control (WT) : It is an important reference mark for whether the PCR reaction is successful and whether the product band position and size meet the theoretical requirements.

② Control (B) : PCR amplification was performed without template in the PCR reagent to monitor whether the reagent was contaminated.

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component	Volume (μl)	
1	2 × Rapid Taq Master Mix (Vazyme P222)	12.5	
2	ddH ₂ O	9.5	
3	Primer A(10pmol/μl)	1	
4	Primer B(10pmol/μl)	1	
5	Template(≈100ng/μl)	1	

PCR program ① priority selection

Seg.	Temp.	Time	Cycle
1	95 °C	5min	
2	98 °C	30s	20×
3	65 °C* (-0.5 °C/cycle)	30s	
4	72 °C	45s*	
5	98 °C	30s	20×
6	55 °C*	30s	
7	72 °C	45s*	
8	72 °C	5min	
9	10 °C	hold	

PCR program ② the second choice

Seg.	Temp.	Time	Cycle
------	-------	------	-------



江苏集萃药康生物科技股份有限公司
GemPharmatech Co.,Ltd

1	95°C	5min	
2	98°C	30s	35×
3	58°C*	30s	
4	72°C	45s*	
5	72°C	5min	
6	10°C	hold	

Note*: Annealing temperature and extension time can be determined according to the actual amplification situation and amplification enzyme efficiency.