


## Genotyping Report

|           |            |             |          |                    |          |
|-----------|------------|-------------|----------|--------------------|----------|
| Strain ID | T033687    | Strain Type | KO(Cas9) | Genetic Background | C57BL/6J |
| Designer  | Shiyu Zhou | Gene Name   | Intu     |                    |          |

### 1. Strategy of Genotyping

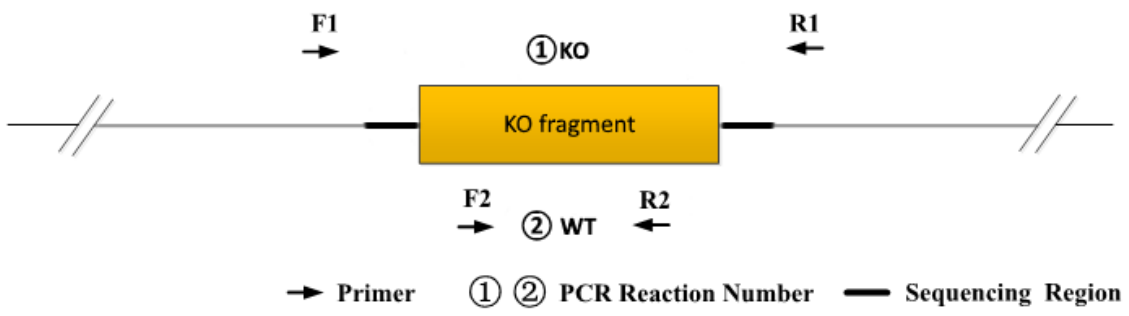


Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.  
Heterozygote: ①PCR reaction obtains a WT band and a KO band; ②PCR reaction obtains a WT band.  
Homozygote: ①PCR reaction obtains a single KO band; ② PCR reaction without product.  
Note: 1)The sizes of WT and Targeted band are shown below.  
2)If the WT band is too large, it may not be possible to obtain a WT band.

### 2. Primer Information

| PCR No. | Primer No. | Primer Name          | Sequence                     | Band Size                      |
|---------|------------|----------------------|------------------------------|--------------------------------|
| PCR①    | F1         | JS24621-Intu-5wt-tF1 | GGTAGACTGGACAGTAGTGTGTGGCTT  | WT:47688 bp<br>Targeted~251 bp |
|         | R1         | JS24621-Intu-3wt-tR1 | CCTGATTGCCCAACACAACCTCCAT    |                                |
| PCR②    | F2         | JS34621-Intu-wt-F1   | G TTCCTATTGTGGGGAGGATAGTCCTG | WT:402 bp<br>Targeted:0bp      |
|         | R2         | JS34621-Intu-wt-R1   | GGAAGCCAAGATGCATAATAGAGCC    |                                |

### 3. Gel Image

|                                                                                          |
|------------------------------------------------------------------------------------------|
| ccagcactggtgacgcttgaagccttagagccct-----47420bp-----AAGaggattcctacttgcttccccactattgatgtct |
|------------------------------------------------------------------------------------------|

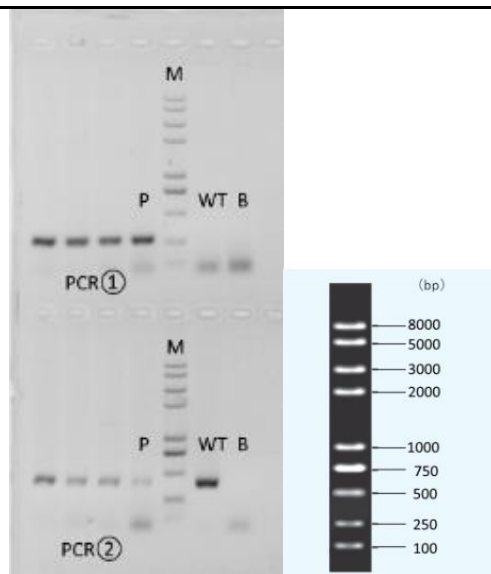

Note: P: Positive control; WT: Wildtype control; B: Blank control (ddH<sub>2</sub>O); M: DNA Ladder

#### 4. PCR Condition

| PCR Reaction Component |                                                 |      |             |
|------------------------|-------------------------------------------------|------|-------------|
| Seg.                   | reaction component                              |      | Volume (μl) |
| 1                      | 2 × Taq Master Mix , Dye Plus, (Vazyme P112-03) |      | 12.5        |
| 2                      | ddH2O                                           |      | 9.5         |
| 3                      | Primer A(10pmol/μl)                             |      | 1           |
| 4                      | Primer B(10pmol/μl)                             |      | 1           |
| 5                      | Template(≈100ng/μl)                             |      | 1           |
| PCR program            |                                                 |      |             |
| Seg.                   | Temp.                                           | Time | Cycle       |
| 1                      | 95℃                                             | 5min | 20×         |
| 2                      | 98℃                                             | 30s  |             |
| 3                      | 65℃ （-0.5℃/cycle）                               | 30s  |             |
| 4                      | 72℃                                             | 45s  |             |
| 5                      | 98℃                                             | 30s  | 20×         |
| 6                      | 55℃                                             | 30s  |             |
| 7                      | 72℃                                             | 45s  |             |
| 8                      | 72℃                                             | 5min |             |
| 9                      | 10℃                                             | hold |             |

Reviewer:Ting Sun

Date:2020-2-26